(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization International Bureau



(43) International Publication Date 29 August 2002 (29.08.2002)

PCT

(10) International Publication Number WO 02/066973 A2

(51) International Patent Classification⁷: G01D 18/00 G01N 33/00,

(21) International Application Number: PCT/US02/01294

(22) International Filing Date: 15 January 2002 (15.01.2002)

(25) Filing Language:

English

(26) Publication Language:

English

(30) Priority Data: 09/784,710 15 February 2001 (15.02.2001) US

(71) Applicant: ABBOTT LABORATORIES [US/US]; D-377 AP6D, 100 Abbott Park Road, Abbott Park, IL 60064-6050 (US). A.; 2467 Higgins Place, Dublin, OH 43016 (US). BORUSZEWSKI, Dennis, V.; 542 Warwick Lane, Pickerington, OH 43147 (US). GUNTRIP, David, P.; 7326 Cook Farm Road, Reynoldsburg, OH 43068 (US).

(74) Agents: WEINSTEIN, David, L. et al.; D-377 AP6D, 100 Abbott Park Road, Abbott Park, IL 60064-6050 (US).

(81) Designated States (national): CA, JP, MX, NO.

(84) Designated States (regional): European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR).


Published:

 without international search report and to be republished upon receipt of that report

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

(72) Inventors: WANG, Jianjun; 8241 Glencullen Court, Dublin, OH 43017 (US). MONDIEK, David,

(54) Title: METHOD AND APPARATUS FOR CALIBRATION OF INSTRUMENTS THAT MONITOR THE CONCENTRATION OF A STERILANT IN A SYSTEM

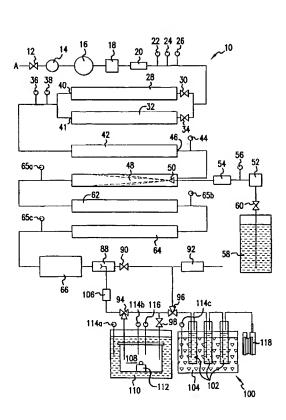

(57) Abstract: A method and apparatus for calibrating a sensor for determination of the concentration of a sterilant, e. g., hydrogen peroxide vapor, in a sterilization system. This invention provides a method for calibrating a sensor that is used for measuring the quantity of a sterilant in a system for delivering the sterilant, the method comprising the steps of: (a) generating reference calibration data, the reference calibration data showing a mathematical relationship between a measurable parameter, e. g. voltage, and a quantity of the sterilant, e. g., parts of sterilant per million parts of air (ppm), for a plurality of sensors; (b) generating sensor calibration data, the sensor calibration data showing a mathematical relationship between the measurable parameter and the quantity of the sterilant for an individual sensor; and(c) normalizing the sensor calibration data to compensate for the difference between the measurable parameter for the reference calibration data and the measurable parameter for the sensor calibration data, whereby data obtained by the individual sensor can be used to accurately determine the quantity of sterilant in the system.

WO 02/066973 A2

METHOD AND APPARATUS FOR CALIBRATION OF INSTRUMENTS THAT MONITOR THE CONCENTRATION OF A STERILANT IN A SYSTEM

BACKGROUND OF THE INVENTION

5

10

1. Field of the Invention

This invention relates to a method for calibration of instruments. More particularly, this invention relates to a method for calibration of instruments that monitor the concentration of a sterilant, e. g., hydrogen peroxide, in a system.

2. Discussion of the Art

15

20

25

Aseptic processing of consumable products, such as nutritional compounds and food products, is typically effected by separate sterilization of the products and the containers within which the products are packaged. Subsequent to sterilization, the sterilized products are placed in sterilized containers and sealed in a sterile environment for shipment, storage, and use.

Sterilization of such containers, which may include sterilization of separate closures as well, can be performed efficiently by use of a sterilant, such as hydrogen peroxide (H_2O_2) vapor, prior to the introduction of the desired sterilized products into the containers. In such a process, the containers are introduced into a sterilization apparatus in which the containers are flushed with hydrogen peroxide vapor. The containers are subsequently flushed with warm air or any other fluid suitable for achieving desirably low levels of residual hydrogen peroxide. This general procedure is highly effective in achieving sterilization of the containers and can be performed on any other articles that will come into contact with the material to be introduced into the containers.

30 containe

Notwithstanding the effectiveness of sterilization by hydrogen peroxide, accurate monitoring of concentration levels of hydrogen peroxide vapor can be problematic. Problems in monitoring the concentration of hydrogen peroxide vapor result in part from changes in the physical and chemical properties of hydrogen

peroxide vapor under processing conditions and the decomposition of hydrogen peroxide vapor upon contact with surfaces of various objects within the processing area. As such, undesired deviations of the concentration of hydrogen peroxide vapor from a process set point, along with excessive decomposition of hydrogen peroxide vapor, can result in loss of sterility of the containers and surrounding aseptic processing area. Moreover, hydrogen peroxide vapor is corrosive in nature, and thus excessive concentration levels of hydrogen peroxide may bring about detrimental effects to the equipment in and surrounding the processing area and surfaces of objects within the processing area. Furthermore, in accordance with government standards, subsequent use of sterilized containers requires low levels of residual sterilant.

Heretofore, detection systems for hydrogen peroxide vapor have been undesirably bulky, as exemplified by conventional near infrared (NIR) analysis apparatus. In addition, current off-line testing methods are typically too slow for monitoring levels of sterilant with sufficient accuracy. Previous arrangements have not allowed real time monitoring throughout an aseptic processing cycle, and in particular, have not been capable of monitoring concentrations of hydrogen peroxide vapor within the sterilization apparatus at select locations along the sterilant supply system during actual operations. However, U. S. Patent No. 5,608,156 and Taizo et al., "Application of a Newly Developed Hydrogen Peroxide Vapor Phase Sensor to HPV Sterilizer", PDA Journal Of Pharmaceutical Science & Technology, Vol. 52, No. 1/January-February 1998, pp. 13-18, describe methods of detecting the concentration of hydrogen peroxide vapor and an apparatus therefor that appear to address some of the foregoing problems.

The concentration of sterilants detected within a system is generally a function of various environmental parameters, such as, for example, temperature, relative humidity, and various measurement conditions, such as, for example, proximate location of measurement. Conventional detection systems for sterilant typically cannot or do not account for fluctuations of environmental parameters and measurement conditions. However, such fluctuations can substantially affect the results of signal generation and data collection when commercially available sensors and equipment are used. It is therefore beneficial to maintain operating parameters proximate the location of measurement as uniform as possible during data collection.

United States Serial No. 09/443,768, filed November 9, 1999, entitled STERILANT MONITORING ASSEMBLY AND APPARATUS AND METHOD USING SAME, incorporated herein by reference, describes an integrated system for determining the concentration of hydrogen peroxide for aseptic process validation, control, and monitoring. This system is compact and can be used for on-line determination of the concentration of hydrogen peroxide. The system requires a unique calibration procedure at regular intervals to guarantee reliable and accurate test results. This system utilizes a sensor having elements made of SnO₂. When SnO₂ is heated to a high temperature, around 400° C, in the absence of oxygen, free electrons flow easily through the grain boundary of the SnO₂ particles. In clean air, oxygen, which traps free electrons by its electron affinity, is adsorbed onto the surface of the SnO₂ particle, forming a potential barrier in the grain boundaries that restricts the flow of electrons, thereby causing the electronic resistance to increase. When the sensor is exposed to hydrogen peroxide vapor, SnO₂ adsorbs its gas molecules and causes oxidation. This lowers the potential barrier, allowing electrons to flow more easily, thereby reducing the electrical resistance. Thus, the sensor uses an indirect method to measure the concentration of hydrogen peroxide vapor.

10

15

20

25

30

Voltage data from the output of the sensor must be compared to a database derived from a calibration process. The output of two different sensors cannot be compared directly without calibration. The calibration procedure uses several representative points (i. e., concentration at a given voltage) to establish a mathematical relationship that covers a specific test window. Only by means of calibration can the output voltage of a sensor be converted to a value of concentration.

Calibration procedures are important for minimizing deviations caused by such components as semiconductor chips, batteries, and signal conditioning circuits in a sensor in a portable detection system. Calibration procedures are important for minimizing deviations caused by such components as temperature and humidity compensation circuits, heating coils, data recording systems, and memory chips in a sensor in a fixed detection system.

If the calibration method is not reliable, the concentration of hydrogen peroxide vapor detected by a sensor might be misleading. In turn, an erroneous determination of the concentration of hydrogen peroxide vapor can bring about contamination in the operation system and result in spoilage. For example, a drop in voltage in the

response of the sensor caused by an increase in the rate of flow of air may be interpreted as a decrease in the concentration of hydrogen peroxide vapor in the system. This apparent decrease may cause the controls in the system to increase the quantity of hydrogen peroxide delivered, thereby providing an excessive amount of hydrogen peroxide vapor. An excessive amount of hydrogen peroxide in the system may result in an excessive amount of residue. Conversely, an increase in voltage in the response of the sensor may result from a decrease in the rate of flow of air. If the delivery rate of hydrogen peroxide is correspondingly reduced, a breach in the sterility of the system may occur.

Calibration of sensors one at a time is inefficient, and, consequently, costly. It is well-known that no two sensors chosen at random are likely to be identical.

Accordingly, it would be desirable to find way to calibrate individual sensors accurately and at reasonable cost. In addition, it would be desirable to find a way to calibrate individual sensors so that one or more of them could be used in portable units. The use of a greater number of portable units is desirable so that measurement of the concentration of hydrogen peroxide can be made at any point in a production line.

SUMMARY OF THE INVENTION

20

25

30

5

10

15

This invention provides a method and apparatus for calibrating a sensor for determination of the concentration of a sterilant, e. g., hydrogen peroxide vapor, in a sterilization system.

In one aspect, this invention provides a method for calibrating a sensor that is used for measuring the quantity of a sterilant in a system for delivering the sterilant, the method comprising the steps of:

- (a) generating reference calibration data, the reference calibration data showing a mathematical relationship between a measurable parameter, e. g. voltage, and a quantity of the sterilant, e. g., parts of sterilant per million parts of air (ppm), for a plurality of sensors;
- (b) generating sensor calibration data, the sensor calibration data showing a mathematical relationship between the measurable parameter and the quantity of the sterilant for an individual sensor; and

(c) normalizing the sensor calibration data to compensate for the difference between the measurable parameter for the reference calibration data and the measurable parameter for the sensor calibration data, whereby data obtained by the individual sensor can be used to accurately determine the quantity of sterilant in the system.

The reference calibration data can be generated by a method comprising the steps of:

(a) providing a plurality of sensors;

10

15

20

5

- (b) subjecting each of the plurality of sensors to at least two quantities of air (e. g., 30 cubic meters/hour and 110 cubic meters/hour), each of the at least two quantities of air having (1) a known quality (e. g., 60% relative humidity at 70 °C) and (2) a known concentration of sterilant vapor (e. g., 10,000 ppm of hydrogen peroxide vapor), the sterilant vapor having a known physical condition (e. g., 70 °C);
- (c) measuring the signals (e. g., voltage) emitted by each of the plurality of sensors, each of the signals being proportional to a concentration (e. g., ppm) of sterilant vapor;
- (d) establishing a mathematical relationship between the signals emitted by each of the plurality of sensors and the concentrations of sterilant vapor for each of the plurality of sensors; and
- (e) establishing the reference calibration data by means of a statistical analysis of the signals emitted by each of the plurality of sensors and the concentrations of sterilant vapor for each of the plurality of sensors.

25

The sensor calibration data can be generated by a method comprising the steps of:

(a) providing a sensor;

30

(b) subjecting the sensor to at least two quantities of air, each of the at least two quantities of air having (1) a known quality and (2) a known concentration of sterilant vapor, the sterilant vapor having a known physical condition;

(c) measuring the signals (e. g., voltage) emitted by the sensor, each of the signals corresponding to a concentration of sterilant vapor (e. g., ppm of hydrogen peroxide vapor); and

(d) establishing a mathematical relationship between the signals emitted and the concentrations of sterilant vapor for the sensor.

The sensor calibration data can be normalized to compensate for the difference between the measurable parameter for the reference calibration data (e. g., voltage) and the measurable parameter for said sensor calibration data for the sensor (e. g., voltage) by a method comprising the steps of:

(a) selecting a concentration of sterilant vapor;

5

10

15

20

25

30

- (b) determining the value of the measurable parameter at which the concentration of the sterilant vapor obtained from the sensor calibration data equals the concentration of said sterilant vapor obtained from the reference calibration data; and
- (c) adjusting the values measured by the individual sensor a sufficient amount to compensate for the deviation between the reference calibration data and the sensor calibration data.

The method of this invention brings about a reduction in the time required to calibrate a sterilization system and provides a means for directly comparing signals obtained from different sensors so that the complex steps involved in the indirect measurement of parameters are avoided.

In one embodiment, which is relatively simple to implement, a linear relationship between the measured signal (e. g., voltage) and the concentration of sterilant (e. g., ppm of hydrogen peroxide vapor) is assumed. Points on the curve for an individual sensor (i. e., sensor calibration data) are moved either vertically or rotationally or both in order to convert the values of concentration associated with that point to a value of concentration on the curve representing the reference calibration data.

The form of the data for calibration of sensors for determining the concentration of hydrogen peroxide vapor in a stream of air is not critical. Preferred forms for

presentation of data include, but are not limited to, curves plotted on Cartesian coordinates, nomograms, and look-up tables containing signal/concentration data.

In one embodiment, in which curves plotted on Cartesian coordinates are employed, the method involves the steps of

5

- (a) preparing a reference calibration curve, the reference calibration curve having a slope and an intercept;
- (b) preparing a sensor calibration curve, the sensor calibration curve having a slope and an intercept; and

10

15

(c) normalizing the sensor calibration curve to compensate for (1) the difference between the slope of the reference calibration curve and the slope of the sensor calibration curve and (2) the difference between the intercept of the reference calibration curve and the intercept of the sensor calibration curve.

For this embodiment, the reference calibration curve is prepared by a method comprising the steps of:

(a) providing a plurality of sensors;

20

(b) subjecting each sensor of the plurality of sensors to at least two quantities of air, each of the at least two quantities of air having (i) a known quality and (ii) a known concentration of sterilant vapor, the sterilant vapor having a known physical condition;

(c) measuring the signals emitted by each of the plurality of sensors, each of the signals being proportional to the concentration of sterilant vapor;

25

(d) establishing a linear mathematical relationship between the signals emitted by each of the plurality of sensors and the concentrations of sterilant vapor for each of the plurality of sensors; and

(e) establishing the reference calibration curve by means of a statistical analysis of the signals emitted by each of the plurality of sensors and the concentrations of sterilant vapor for each of the plurality of sensors.

30

For this embodiment, the sensor calibration curve is prepared by a method comprising the steps of:

(a) providing an individual sensor;

5

10

15

20

25

30

- (b) subjecting the sensor to at least two quantities of air, each of the at least two quantities of air having (1) a known quality and (2) a known concentration of sterilant vapor, the sterilant vapor having a known physical condition;
- (c) measuring the signals emitted by the sensor, each of the signals being proportional to the concentration of sterilant vapor; and
- (d) establishing a linear mathematical relationship between the signals emitted by the individual sensor and the concentrations of sterilant vapor for the individual sensor.

The sensor calibration curve is normalized to compensate for (1) the difference between the slope of the reference calibration curve and the slope of the sensor calibration curve and (2) the difference between the intercept of the reference calibration curve and the intercept of the sensor calibration curve by a method comprising the steps of:

- (a) determining the intercept of the reference calibration curve;
- (b) determining the intercept of the sensor calibration curve;
- (c) determining the slope of the reference calibration curve;
- (d) determining the slope of the sensor calibration curve;
- (e) adjusting the sensor calibration curve, if necessary, in order to compensate for the difference between the intercept of the reference calibration curve and the intercept of the sensor calibration curve; and
- (g) adjusting the sensor calibration curve, if necessary, in order to compensate for the difference between the slope of the reference calibration curve and the slope of the sensor calibration curve.

In another aspect, this invention provides a method for calibrating a portable unit for measuring the concentration of hydrogen peroxide. In this method, hydrogen peroxide vapor can be passed through a calibration vessel, which is submerged in a water bath, under controlled testing conditions. A portable sensor for the detection of hydrogen peroxide vapor installed within the calibration vessel will respond to the concentration of hydrogen peroxide, the temperature, and the relative humidity for

each test run. The concentration of hydrogen peroxide can be rechecked by a standard titration method. The accuracy of the value of the concentration of hydrogen peroxide thus determined can be checked by means of a detection unit for hydrogen peroxide vapor residue (Drager kit) to determine the concentration of residual hydrogen peroxide in the flow stream exiting the system. Additional impingers can be added to the system to capture all hydrogen peroxide from the stream of flowing air if residual hydrogen peroxide should be detected by a detection unit for hydrogen peroxide vapor residue. Additional data can be generated to cover those conditions that would be expected to be encountered in a practical application.

5

10

15

20

25

30

The method of this invention is capable of conducting calibrations that mimic both static and dynamic processing conditions, such as seen in sterilization of Bosch aseptic machines and hydrogen peroxide spray sterilization of bottles. Flexibility in the selection of parameters encountered in the generation of hydrogen peroxide vapor allows dynamic and static calibration of instruments for the determination of the concentration of hydrogen peroxide under conditions of high temperature, high concentration, variable rate of flow of air, and variable humidity for the test cell.

The method of this invention allows the calibration of instruments for determining the concentration of hydrogen peroxide vapor for numerous applications. The combined application of a water bath, a titration station, a Drager kit, a metering pump, and a hydrogen peroxide flow meter allows calibration of a portable detection unit for the determination of low concentrations of hydrogen peroxide.

The method of this invention can also be used to investigate the condensation of hydrogen peroxide by changing the rate of injection of hydrogen peroxide, the temperature(s) of the heat exchanger(s), the temperature(s) of the water bath(s), and the rates of flow of air. The rates of condensation and decomposition of hydrogen peroxide can be calculated by determining the difference between the amount of hydrogen peroxide injected and the amount of hydrogen peroxide detected over a selected period of time. The temperature at which air can be saturated by hydrogen peroxide under known conditions can be estimated as well.

The method of this invention can be used to calibrate a sensor for a sterilization system to improve the ability to determine the concentration of sterilant at a given point in the system at any time at various processing stages. By this method, accurate information on the process at each stage can be provided.

The method of this invention can also be used to calibrate instruments for monitoring the concentration of hydrogen peroxide vapor in any environment where hydrogen peroxide vapor is employed. Such environments include, but are not limited to, clean room operations, pharmaceutical isolator sterilization, aseptic processing systems, and microbiological investigation of death rate for various bacteria under the effect of hydrogen peroxide.

5

10

15

25

30

This invention provides a method for calibrating instruments for determining the concentration of a sterilant, e. g., hydrogen peroxide vapor, for a wide range of processing and operating conditions. These operating condition include, but are not limited to, air temperature, flow rate of air, evaporation temperature, hydrogen peroxide injection rate, degree of saturation of air by hydrogen peroxide vapor at different temperatures, humidity variations, pressure fluctuations, etc.

It is to be understood that the method of this invention is not limited to the sterilant hydrogen peroxide. The method of this invention can be used to calibrate instruments for determining the concentration of other gases used for sterilization.

BRIEF DESCRIPTION OF THE DRAWINGS

- FIG. 1 is a schematic diagram that illustrates a system suitable for generating hydrogen peroxide vapor under controllable operating conditions. This diagram includes the components of the system that are needed to calibrate a sensor for a sterilization system.
 - FIG. 2 is a schematic diagram that illustrates a test cell for use in calibration of a sensor for a sterilization system.
 - FIG. 3 illustrates a type of graph suitable for calibrating a sensor for use for detecting the concentration of hydrogen peroxide vapor in a system. The graph includes a reference calibration curve and a sensor calibration curve.

DETAILED DESCRIPTION

5

10

15

20

25

30

As used herein, the term "sterilant" means a substance that can kill one or more target microorganisms. The expression "reference calibration data" means set of data, derived from examining a plurality of sensors, used to construct a database. The set of data is obtained from at least the minimum number of a plurality of sensors that will provide a statistically significant result. The expression "measurable parameter" means a variable that can be measured by a method of measurement acceptable to one of ordinary skill in the art. The term "quantity" means the property or aspect of a thing that can be measured, counted, or compared, such as, for example, rate of flow of air. The expression "sensor calibration data" means set of data, derived from examining an individual sensor, used to construct a database. The set of data is obtained from a single sensor at a given physical condition and at a given concentration level. The term "normalizing" means causing the sensor calibration data to conform to the reference calibration data. The expression "air quality" means the physical condition of the air into which hydrogen peroxide is vaporized. Parameters that can be used to characterize air quality include, but are not limited to, temperature, relative humidity, pressure, and the like. The term "signal" means a fluctuating electric quantity, as voltage, current, electric field strength, whose variations represent coded information. The expression "mathematical relationship" means a mathematical equation that expresses at least one parameter as a function of at least one other parameter. The term "slope" means the rate at which an ordinate of a point of a line on a coordinate plane changes with respect to a change in the abscissa. The term "intercept" means the distance from the origin of coordinates along a coordinate axis to the point at which a line, a curve, or a surface intersects the axis. The expression "Drager kit" means a unit for detecting the concentration of residual hydrogen peroxide vapor. The expression "static condition" refers to the situation in which both the sensor and the hydrogen peroxide vapor are momentarily sealed in an air-tight chamber maintained at a constant temperature and humidity during a test run or calibration run. The expression "dynamic condition" refers to the situation in which the sensor is installed in the stream in which the hydrogen peroxide vapor flows.

In one aspect, this invention involves a method for calibrating a sensor that is used for measuring the quantity of a sterilant in a system for delivering the sterilant, the method comprising the steps of:

5

10

15

20

25

30

(a) generating reference calibration data, the reference calibration data showing a relationship between a measurable parameter, e. g. voltage, and a quantity of the sterilant, e. g., parts of sterilant per million parts of air (ppm), for a plurality of sensors;

- (b) generating sensor calibration data, the sensor calibration data showing a relationship between the measurable parameter and the a quantity of the sterilant for an individual sensor;
- (c) normalizing the sensor calibration data to compensate for the difference between the measurable parameter for the reference calibration data and the measurable parameter for the sensor calibration data for the individual sensor, whereby the quantity of sterilant in the system can be accurately determined.

In order to generate the reference calibration data, each of a plurality of sensors is subjected to at least two quantities of air (e. g., 30 cubic meters/hour and 110 cubic meters/hour). Each of the at least two quantities of air has (1) a known quality (e. g., 60% relative humidity at 70 °C) and (2) a known concentration of sterilant vapor (e. g., 10,000 ppm of hydrogen peroxide vapor). The sterilant vapor also has a known physical condition (e.g., 70 °C). The signals emitted by each of the sensors are measured. Each of the signals measured corresponds to the concentration of sterilant vapor. Then, a mathematical relationship is established between the signals measured and the concentrations corresponding to the measured signals for each of the sensors. The foregoing operations are repeated until a mathematical relationship between the signals measured and the concentrations corresponding to the measured signals has been established for each sensor of the plurality of sensors. Then, the reference calibration data is derived from the data and mathematical relationships obtained from the foregoing measurements. Typically, the reference calibration data represents the mean (average) of the data obtained from the foregoing measurements. The number of sensors used to prepare the reference calibration data is preferably sufficient to provide reference calibration data that exhibits a normal distribution with about 95% of the values of the sensors falling within two standard deviations of the mean. In general, it is preferred that at least 10 sensors be used to provide data for preparing the reference calibration data.

In order to generate the sensor calibration data, an individual sensor is provided. This sensor is subjected to at least two quantities of air, each of the at least two quantities of air having (1) a known quality and (2) a known concentration of sterilant vapor. The sterilant vapor has a known physical condition. The signals emitted by the sensor are measured. Each of the measured signals corresponds to a concentration of sterilant vapor. A mathematical relationship is established between the signals emitted by the sensor and the concentrations corresponding to the signals measured for the sensor. Then, the sensor calibration data is derived from the data and mathematical relationships obtained from the foregoing measurements

Then the sensor calibration data is normalized to compensate for the difference between the measurable parameter for the reference calibration data (e. g., voltage) and the measurable parameter for the sensor calibration data for the individual sensor. In order to carry out the normalization step, a concentration of sterilant vapor is selected. The value of the measurable parameter at which the concentration of the sterilant vapor obtained from the sensor calibration data equals the concentration of said sterilant vapor from the reference calibration data is determined. This value of the measurable parameter for a given concentration must be in both the reference calibration data and the sensor calibration data. The values of the measurable parameter in the sensor calibration data are adjusted a sufficient amount to compensate for the deviation of the sensor calibration data from the reference calibration data.

10

15

20

25

30

The form of the data for both the reference calibration data and the sensor calibration data is not critical. Preferred forms of data include, but are not limited to, curves plotted on Cartesian coordinates, nomograms, and tables containing signal/concentration data. Of course, the form for the reference calibration data and the sensor calibration data must be the same in order to ensure a meaningful adjustment.

In one embodiment, in which curves plotted on Cartesian coordinates are employed, the calibration can be carried out in the following manner. A linear reference calibration curve having a slope and an intercept is prepared. The value of the concentration is plotted on the abscissa and the value of signal is plotted on the ordinate. A linear sensor calibration curve for an individual sensor having a slope and an intercept is prepared. The value of the concentration is plotted on the abscissa and the value of signal is plotted on the ordinate. The reference calibration curve and

the sensor calibration curve intersect at a set of coordinates. The sensor calibration curve is normalized to the reference calibration curve to compensate for (1) the difference between the slope of the reference calibration curve and the slope of the sensor calibration curve and (2) the difference between the intercept of the reference calibration curve and the intercept of the sensor calibration curve. Normalization can be carried out in several ways. In one embodiment of normalization, the sensor calibration curve can be adjusted physically, by adjusting it vertically and rotationally so that the intercept and the slope of the sensor calibration curve conform to the intercept and the slope, respectively, of the reference calibration curve. In another embodiment of normalization, the variable representing the signal on the sensor calibration curve (i. e., voltage or log (voltage)) can be substituted into the equation that relates the concentration and the signal on the reference calibration curve.

5

10

15

20

25

30

Referring now to FIGS. 1 and 2, compressed air at room temperature is introduced into a system 10 suitable for generating hydrogen peroxide vapor under controllable operating conditions by means of a valve 12. The compressed air is represented by the letter "A". The air is filtered by means of an air filter 14 to remove dirt and other foreign substances. A valve 16, preferably a diaphragm valve, is used to control the flow rate of the air, and a pressure regulator 18 is used to adjust the pressure of the air. The rate of flow of the air is monitored by a flow meter 20 and the parameters of air quality are measured and recorded by a data acquisition system (e. g., Keithley 500 data acquisition system, Keithley Instruments, Inc., 28775 Aurora Rd., Cleveland, Ohio 44139). The parameters of air quality include temperature, humidity, and pressure. These parameters are monitored by a temperature sensor 22, a humidity sensor 24, and a pressure sensor 26, respectively. Because of seasonal variations in the quality of air in the room containing the components of the system 10, the moisture content of the air is adjusted to a predetermined level for the purpose of the calibration. The moisture content of the air is adjusted by using either a dehumidifier 28, which is operated by opening or closing a valve 30, or a humidifier 32, which is operated by opening or closing a valve 34. The adjusted air quality is measured by means of a temperature sensor 36 and a humidity sensor 38 located downstream of the exit port 40 of the dehumidifier 28 and the exit port 41 of the humidifier 32. Accordingly, the temperature level and the humidity level of the air will be known.

The air is then heated to a specified temperature by means of a heat exchanger 42. The concentration of hydrogen peroxide (hereinafter H₂O₂) vapor in a sample of air saturated with H₂O₂ vapor is a function of temperature. As the temperature of the air is increased, the air can hold a higher concentration of H₂O₂ vapor. The temperature of the air can be monitored by a temperature sensor 44 located downstream of the exit 46 of the heat exchanger 42. In the heat exchanger 42, temperature is a function of steam pressure. By adjusting the steam pressure of the heat exchanger 42, the temperature of the air can be controlled. The air is then introduced into a heat exchanger 48 and merged with a stream of H₂O₂ droplets created by an atomizer 50. The H₂O₂ is caused to evaporate rapidly and is further heated to a higher temperature. A metering pump 52 controls the rate of injection of liquid H_2O_2 into the air flowing into the heat exchanger 48. A flow meter 54 is used to monitor the flow rate of liquid H₂O₂ to ensure that the rate of injection of H₂O₂ is correct. The temperature of the liquid H₂O₂ can be measured by a temperature sensor 56 to ensure a consistent temperature for all calibration operations. A tank for storing the liquid H₂O₂ is designated by the reference numeral 58, and a valve for controlling the injection of liquid H₂O₂ into the stream of air is designated by the reference numeral 60. The tank 58 and the valve 60 are those conventionally used for introducing liquid H₂O₂ into a sterilization system. Accordingly, the tank 58 and the valve 60 are resistant to H₂O₂.

10

15

20

25

30

The mixture of air and H_2O_2 vapor may be in one of several conditions. These conditions may be described as follows:

- 1. Air is saturated with H_2O_2 vapor condensation of droplets of the H_2O_2 vapor occurs when temperature decreases;
- 2. Air is not saturated with H_2O_2 vapor H_2O_2 vapor decomposes on the surfaces of certain materials; H_2O_2 vapor condenses when temperature drops below the saturation temperature;
- 3. Under both conditions (saturated and unsaturated), if the temperature increases, H₂O₂ droplets will revaporize;
- 4. Latent heat of evaporation of H₂O₂ is positive; evaporation of H₂O₂ generates extra heat in the system.

The temperature of the air is generally higher than the temperature at which the heat exchanger 48 is set, and the temperature of the air is difficult to predict and control.

Two additional heat exchangers 62 and 64 can be used to reduce the temperature of the air and ensure that the air is saturated with H_2O_2 vapor at a known temperature. For the air to reach a state of saturation, it is common to evaporate a liquid at a high temperature and then reduce the temperature. If a sufficient quantity of liquid is injected into the system at a high temperature, saturation can be guaranteed when the vapor is cooled. Also, the saturation process should be conducted at the last possible point before the H_2O_2 vapor is to be used for sterilization, because it is difficult to control condensation and decomposition at some other point in the system.

5

10

15

20

25

30

The temperature of the air at various locations in the system 10 can be monitored by temperature sensors 65a, 65b, and 65c.

The air containing the H_2O_2 is then introduced into a test cell 66. A simplified diagram of the test cell 66 is shown in FIG. 2. The test cell 66 is described in U. S. Serial No. 09/443,768, filed November 19, 1999, entitled STERILANT MONITORING ASSEMBLY AND APPARATUS AND METHOD USING SAME, incorporated herein by reference. In U. S. Serial No. 09/443,768, the test cell is referred to as a sterilant monitoring assembly and is designated by the reference numeral 200. A sensor 68 suitable for determining the concentration of H_2O_2 vapor is installed in the test cell 66 in order to be calibrated. It is preferred that the conditions within the test cell 66 be substantially identical to those used in an actual sterilization process. It should be noted that the data for the calibration operation is likely to be valid only after the system 10 reaches a state of equilibrium. Accordingly, the system 10 should be run a sufficient amount of time in order to reach a state of equilibrium.

It is preferred that two types of runs be performed in each calibration operation. In a static test run, an internal cell 70 in the test cell 66 is placed in fluidic communication with the flowing stream of air for a specified amount of time, the flowing stream of air having reached a state of equilibrium. In order to determine the concentration of H_2O_2 vapor in a sample of air, the internal cell 70 in the test cell 66 is then sealed by means of a piston 72 moving in a cylinder 74 under the force of compressed air pushing against a spring 76. The compressed air can be introduced to the cylinder 74 through an inlet 78. An H_2O_2 -resistant O-ring 80 can be used to ensure that the internal cell 70 is sealed against a heat isolation plate 82. The cover

of the test chamber 66 is designated by the reference numeral 84 and a thermocouple is represented by the reference numeral 86. The data obtained by means of this static test run reflects the level of concentration of H_2O_2 vapor, independent of the rate of flow of air. The flowing air surrounding the internal cell 70 serves as an insulating jacket to maintain a uniform and constant temperature during the static test run. In a dynamic test run, the internal cell 70 is placed in fluidic communication with the flowing air at all times. In other words, the internal cell 70 is not sealed against the heat isolation plate 82. In the dynamic test run, the data obtained only reflect the concentration of H_2O_2 vapor at a specific flow rate of air.

10

15

20

25

30

The results obtained from the two types of calibration runs can be used to verify the determination of the concentration of H₂O₂ vapor under actual sterilization conditions. For example, if the readings of concentration of H₂O₂ vapor determined in static test runs remain substantially constant, but if the readings of concentration of H₂O₂ vapor determined in dynamic test runs vary substantially, (1) the flow rate of the air may have changed, (2) the rate of addition of sterilant to the system may have changed, or (2) something may be wrong with the system. Similarly, if the readings of concentration of H₂O₂ vapor determined in static runs conditions vary substantially. but if the readings of concentration of H₂O₂ vapor determined in dynamic test runs remain substantially constant, (1) the flow rate of the air may have changed, (2) the rate of addition of sterilant to the system may have changed, or (2) something may be wrong with the system. Assuming that the flow rate of the air has not changed, and that the rate of addition of sterilant to the system has not changed, what may be wrong with the system is that the sensor for determining the concentration of H₂O₂ vapor is providing unreliable readings, which may be the result of unreliable calibration. If the method for calibrating the sensor is unreliable, the concentration of H₂O₂ vapor determined by the sensor may be misleading. An erroneous determination of the concentration of H₂O₂ vapor can result in contamination of the system and result in spoilage of product. For example, a drop in voltage in the response measured by the sensor caused by an increase in the rate of flow of air may be interpreted as a decrease in concentration of H₂O₂ vapor in the system. This apparent decrease in response may cause controls in the system to increase the quantity of H₂O₂ delivered, thereby providing an excessive amount of H₂O₂ vapor. An excessive amount of H₂O₂ in the system will result in an excessive amount of residue. Conversely, an increase in voltage in the response of the sensor may result from a

decrease in the rate of flow of air. If the delivery rate of H_2O_2 is correspondingly reduced, a breach in the sterility of the system may occur. For the foregoing reasons, it is frequently desirable to cross-check calibration runs of sensors.

As stated previously, the system 10 should be run a sufficient amount of time in order to reach a state of equilibrium. When the system is activated initially, the air exiting the test cell 66 flows into a sample collection port 88 positioned downstream from the test cell 66. A valve 90 is operated to allow the air flowing through the sample collection port 88 to be transported to a catalytic converter 92, which accelerates the decomposition of H_2O_2 (i. e., breakdown of H_2O_2) before it enters the ambient environment. During this initial activation period, a valve 94, preferably a programmable three-way valve, allows air to flow through the system 10 to a valve 96, preferably a three-way valve, and from there to the catalytic converter 92. During this initial activation phase, a valve 98, which will be described later, is closed.

10

15

20

25

30

After the initial activation phase, i. e., when the system 10 has reached a state of equilibrium, the signal (e. g., voltage data) obtained from the sensor 68 is correlated with the concentration of H₂O₂ vapor in the test cell 66. In order to perform this correlation, a sample of air collected via the sample collection port 88 is passed through a titration station 100, which comprises a series of impingers 102 set in an ice bath 104. The flow rate of the sample of air is measured by a flow meter 106; the duration of the test is precisely controlled by the operation of the valve 94. In order to perform this measurement, the valve 94 allows air to flow through the system 10 to the valve 96, which is operated to allow the air to enter the titration station 100. The valve 90 is operated to allow air to flow to the catalytic converter 92, and the valve 98 is closed. The amount of H₂O₂ captured in the impingers 102 is a function of time. By measuring the time consumed and the amount of H₂O₂ collected in the impingers 102, the concentration of H_2O_2 can be estimated. For example, if the flow of air is 10 cubic meters per minute, and if approximately 10 mL of H₂O₂ vapor condense in the impingers 102, then the concentration of H₂O₂ is equal to approximately 1 mL per cubic meter. The rate of flow of the sample of air to the titration station 100 can be controlled by adjusting the size of an orifice in the sample collection port 88.

When the titration station 100 is used to determine the concentration of H_2O_2 in the sample of air in the test cell 66, the valve 94 is operated to allow the air containing the H_2O_2 vapor to bypass a test chamber 108 in a water bath 110. The test chamber 108 can be used to calibrate a portable H_2O_2 detection unit 112. A portable H_2O_2

5

10

15

20

25

30

detection unit is one that can be installed in a bottle sterilizer, such as, for example, a Bosch apparatus, or that can be placed at any location in an isolator or a clean room to monitor sterilization by H_2O_2 . A portable H_2O_2 detection unit can also be installed inside a container, such as a 32-ounce plastic container, to detect the H₂O₂ vapor condition after an H₂O₂ vapor spraying procedure for sterilization of bottles. The portable unit is an integrated unit comprising an H₂O₂ sensor, a data recording element, a battery for power supply, a signal conditioning circuit, a near infrared (NIR) receiver, a DC/DC converter, and a temperature sensor. Portable units are preferably programmable, and the data obtained by using such a unit can downloaded after a test. A portable unit has a size of about 1inch by 2 inches by 2 inches. Portable H₂O₂ detection units are described in U. S. Serial No. 09/360,772, filed July 26, 1999, entitled SELF-CONTAINED STERILANT MONITORING ASSEMBLY AND METHOD OF USING SAME, incorporated herein by reference. Before measurements are performed in the test chamber 108, it is preferred that the system 10 be run for a sufficient amount of time to avoid condensation of H₂O₂. Such condensation results from the temperature differential between the air containing the H2O2 vapor and the internal surfaces of pipes, the test cell 66, the test chamber 108, and the valves. During the initial activation phase that is suggested prior to utilizing the test chamber 108, the valve 94 is operated to allow air to flow into the test chamber 108, and the valve 98 is open. The valve 96 is operated to allow the air that exits from the test chamber 108 to flow to the catalytic converter 92, and the valve 90 is open. When the test chamber is used to calibrate a portable detection unit 112, the valve 96 is operated to prevent the air that exits from the test chamber 108 from flowing to the catalytic converter 92, and the valve 90 is open. The valve 96 is operated to allow the air that exits from the test chamber 108 to flow to the titration station 100.

During the calibration of the portable unit 112, the portable unit 112 is placed inside the test chamber 100 and operated by a remote control (not shown), preferably a near infrared (NIR) remote control. The flow of air is adjusted by adjusting the size of the orifice of the sample collection unit 88. The temperature of the water bath 110 is the same as the temperature of the sample of air containing the H_2O_2 vapor. By measuring the duration of the run and the flow rate of the sample of air, the concentration of H_2O_2 vapor in the sample can be calculated in the same manner as was used to calibrate the sensor 68 tested in the test cell 66. A relationship between

the voltage data obtained from the portable H_2O_2 detection unit 112 and the concentration value of H_2O_2 vapor can then be established.

5

10

15

20

25

30

The valve 90 is allowed to remain open most of the time, but is generally closed when a low flow rate of air is required. If the flow rate of air is low, the valve 90 is closed and the air will flow only through the valves 94, 96, and 98. It should be noted that the valves 12, 16, 30, 34, and 60 are operated as required in order to set the conditions of flow rate and humidity and vary the quantity of H_2O_2 admitted to the system 10.

Additional temperatures 114a, 114b, and 114c can be used to monitor temperature in the test chamber 108 and in the titration station 100. A pressure sensor 116 can be used to monitor pressure in the test chamber 108. A Drager kit 118 can be used to detect the concentration of residual H₂O₂ vapor that exits the system 10.

The components of the system including, but limited to air conditioning components (e. g., pressure regulators, flow meters, dehumidifier(s), humidifier(s), heat exchangers), temperature sensors, humidity sensors, atomizer, metering pump, test cell, water bath, test chamber, titration station, and valves, are all commercially available. Optimal selection of these components can be made by one of ordinary skill in the art. While the sizes of the system 10 and the components thereof are not critical, the dimensions of a few components are given to provide an idea of the amount of space required for a calibration system. The size of the system 10 is about 3 meters x 4 meters x 3 meters. The size of the test cell 66 is about 1 foot x 0.5 foot x 0.5 foot. The size of the portable unit 112 is about 1 inch x 2 inches x 2 inches. The size of the Drager kit 118 is about 3 inches x 3 inches x 4 inches. The size of the test chamber 108 is about 8 inches high x about 4 inches in diameter. The volume of an impinger 102 is about 125 mL.

The operating conditions of the system can vary but they generally fall within the following ranges:

- (a) temperature of the air ranges from about -10 °C to about 85 °C;
- (b) relative humidity of the air ranges from about 10 to about 100%;
- (c) the flow rate of air ranges from about 1 to about 200 cubic meters per hour;
 - (d) the concentration of sterilant ranges from 100 to 25000 ppm.

Calibration procedures should be conducted at regular time intervals to obtain reliable and accurate test results. The calibration system 10 can ensure the reliability and accuracy of the prediction of concentration of H_2O_2 determined by an H_2O_2 sensor by simulation of process conditions.

A series of calibration runs can be conducted for different temperature ranges, various H_2O_2 concentration levels, selected flow rates, adjusted relative humidity, different liquid H_2O_2 concentrations.

The method of this invention can be used to calibrate a sensor for a sterilization system to improve the ability to determine the concentration of sterilant at a given point in the system at any time at various processing stages. By this method, accurate information on the process at each stage can be provided.

The method of this invention can also be used to calibrate instruments for monitoring the concentration of H_2O_2 vapor in any environment where H_2O_2 vapor is employed. Such environments include, but are not limited to, clean room operations, pharmaceutical isolator sterilization, aseptic processing systems, and microbiological investigation of death rate for various bacteria under the effect of hydrogen peroxide.

This invention provides a method for calibrating instruments for determining the concentration of a sterilant, e. g., H_2O_2 vapor, for a wide range of processing and operating conditions. These operating condition include, but are not limited to, air temperature, flow rate of air, evaporation temperature, hydrogen peroxide injection rate, degree of saturation of air by hydrogen peroxide vapor at different temperatures, humidity variations, pressure fluctuations, etc.

It is to be understood that the method of this invention is not limited to the sterilant hydrogen peroxide. The method of this invention can be used to calibrate instruments for determining the concentration of other gases used for sterilization.

The following non-limiting example further illustrates the method of this invention.

30

5

10

15

20

25

EXAMPLE 1

The following example illustrates the calibration of a sensor by the method of this invention.

The system 10 shown in FIG. 1 was used to carry out the method of this example. When the system is activated, the valves 12, 16, 30, 34, and 60 are operated so as to deliver the desired quantities of air and H_2O_2 and the desired quality of air to the test cell 66. Because the test cell 66 is being used to calibrate a sensor 66, the valves 94, 96, and 98 are operated so that the air bypasses the test chamber 108. The valve 90 is allowed to remain open during the calibration runs of this example. Typical air quality conditions are 65% relative humidity, 70 °C. Typical flow rates of air range from about 30 cubic meters/hour to about 110 cubic meters/hour.

Reference Calibration Data

Fifteen (15) sensors (TGS 816, Figaro USA, Inc., Glenview, IL 60025) were used to prepare a linear reference calibration curve. This curve is described by the equation $y_1 = a_1 + b_1x$, where x represents concentrations of H_2O_2 , y_1 represents the mean of the voltage readings taken at a particular value of concentration of H_2O_2 for the 15 sensors, a_1 represents the mean (average) of the intercepts of the plots of voltage v. concentration of H_2O_2 for the 15 sensors, and b_1 represents the mean (average) of the slopes of the plots of voltage v. concentration of H_2O_2 for the 15 sensors. See FIG. 3. The value of the slope and the value of the intercept of each of the 15 sensors were determined by data collected at at least three points (concentrations) on the abscissa. Each point was determined by at least two repetitions at each concentration. Table 1 is a compilation of the slopes and intercepts for each of the 15 sensors used to prepare the reference calibration data.

10

15

Table 1

Sensor	Intercept	Slope
1	-0.96	0.39
2	-0.954	0.36
3	-0.938	0.36
4	-1.11	0.41
5	-0.95	0.361
6	-1.04	0.37
7	-1.2	0.34
8	-1.04	0.37
9	-1.3	0.342
10	-1.2	0.45
11	-1.4	0.34
12	-1.11	0.44
13	-0.98	0.361
14	-1.12	0.41
15	-1.05	0.38

The data in Table 1 yielded a mean slope value (i. e., mean of b₁) of 0.379 and a
mean intercept value (i. e., mean of a₁) of -1.090. The following table, Table 2, shows
the relationship between voltage and concentration of H₂O₂ for reference calibration
data utilizing the slope value of 0.379 and the value of -1.090, which were derived
from the data in Table 1.

Table 2

Concentration of H ₂ O ₂ (ppm)	Voltage (volts)	Log (concentration of H ₂ O ₂)	Log (voltage)
3535	1.838	3.548	0.254
5839	2.210	3.766	0.337
8250	2.510	3.916	0.394
10607	2.715	4.026	0.435
12964	2.898	4.113	0.468

Sensor Calibration Data

15

One (1) sensor (TGS 816, Figaro USA, Inc., Glenview, IL 60025) was used to prepare a linear sensor calibration curve. This curve is described by the equation $y_2 = a_2 + b_2x$, where x represents concentrations of H_2O_2 , y_2 represents the voltage reading taken at a particular value of concentration of H_2O_2 for a single sensor, a_2 represents

the intercept of the plot of voltage v. concentration of H_2O_2 for the single sensor, and b_2 represents the slope of the plot of voltage v. concentration of H_2O_2 for the single sensor. See FIG. 3. The values of the slope and the intercept of the single sensor were determined by data collected at at least three points on the abscissa. Each point was determined by at least two repetitions at each concentration. The concentrations used are shown in Table 3.

Table 3

Concentration of H ₂ O ₂ (ppm)	Voltage (volts)	Log (concentration of H ₂ O ₂)	Log (voltage)
3535	1.838	3.548	0.264
5839	2.210	3.766	0.344
8250	2.510	3.916	0.400
10607	2.715	4.026	0.434
12964	2.898	4.113	0.462

The data in Table 3 yielded a slope value of 0.351 and an intercept value of -.098.

Normalizing

The following mathematical relationships were used to convert readings from an individual sensor to standardized values based on the reference calibration data. In the equations that follows, "voltage_r" means voltage of a point on the reference calibration curve, "voltage_s" means voltage of a point on the sensor calibration curve, "concentration _r" means concentration of a point on the reference calibration curve, and concentration _s" means concentration of a point on the sensor calibration curve.

15

5

Reference calibration curve:

$$Log (voltage_r) = -1.090 + 0.379 Log (concentration_r)$$

20 Sensor calibration curve:

$$Log (voltage_s) = -0.981 + 0.352 Log (concentration_s)$$

25

Substituting the equivalent of Log (voltage_s) into the equation for Log (voltage_r) yields:

$$Log (voltage_r) = -1.090 + 0.379 \times [Log (voltage_s] + 0.981]/0.352$$

Thus, one can convert voltage data obtained from an individual sensor to a value of concentration on the reference curve. In order to compare several sensors,

Table 4 shows values encountered in a typical conversion from the value of output voltage of a sensor to a value of voltage on the reference calibration curve, which value corresponds to a specific concentration of H₂O₂ vapor.

Table 4

5

10

15

20

25

Concen- tration of H ₂ O ₂	Log (concen- tration of H ₂ O ₂)	Voltage (ref ¹)	Log (voltage - ref)	Voltage (sen²)	Log (voltage - sen)	Log (voltage converted from voltage (sen))
6000	3.778	2.196	0.342	2.224	0.347	0.342
8000	3.903	2.448	0.389	2.461	0.391	0.389
10000	4.000	2.664	0.426	2.662	0:425	0.426
12000	4.079	2.855	0.456	2.838	0.453	0.456
14000	4.146	3.027	0.481	2.996	0.477	0.481
16000	4.204	3.184	0.503	3.140	0.497	0.503

¹The values of voltage on the reference calibration curve correspond to concentrations obtained by using an equation derived from statistical analysis. The equation was a linear equation in which the slope and an intercept of a line could be determined.

The values of voltage were obtained by an experiment in which the value of voltage corresponded to

The seventh column of Table 4 shows that for a given concentration of H₂O₂, the value of log (voltage) obtained by means of the method of this invention is equal to the value of log (voltage) of the reference calibration curve, as shown in the fourth column of Table 4. Thus, if an individual sensor give a voltage reading of 0.347 volt, the reading will be adjusted to 0.342 volt, which translates to a concentration of H2O2 of 6000 ppm. The method of this invention allows the use of a plurality of sensors to measure the concentration levels of sterilant by means of fixed reference calibration data, thereby reducing the number of data conversions encountered when numerous sensors are used. This feature is extremely useful when several sensors are used at the same time under actual sterilization conditions.

Various modifications and alterations of this invention will become apparent to those skilled in the art without departing from the scope and spirit of this invention, and it should be understood that this invention is not to be unduly limited to the illustrative embodiments set forth herein.

the concentration of H₂O₂.

What is claimed is:

10

15

25

30

A method for calibrating a sensor that is used for measuring the quantity
 of a sterilant in a system for delivering said sterilant, said method comprising the steps
 of:

- (a) generating reference calibration data, said reference calibration data showing a mathematical relationship between a measurable parameter and a quantity of said sterilant, for a plurality of sensors;
- (b) generating sensor calibration data, said sensor calibration data showing a mathematical relationship between said measurable parameter and said quantity of said sterilant for an individual sensor; and
- (c) normalizing said sensor calibration data to compensate for the difference between said measurable parameter for said reference calibration data and said measurable parameter for said sensor calibration data, whereby data obtained by said individual sensor can be used to accurately determine said quantity of sterilant in said system.
- 20 2. The method of claim 1, wherein said reference calibration data is generated by a method comprising the steps of:
 - (a) providing a plurality of sensors;
 - (b) subjecting each of said plurality of sensors to at least two quantities of air, each of said at least two quantities of air having (1) a known quality and (2) a known concentration of sterilant vapor, said sterilant vapor having a known physical condition;
 - (c) measuring signals emitted by each of said plurality of sensors, each of said signals being proportional to a concentration of sterilant vapor;
 - (d) establishing a mathematical relationship between said signals emitted by each of said plurality of sensors and said concentrations of sterilant vapor for each of said plurality of sensors; and
 - (e) establishing said reference calibration data by means of a statistical analysis of said signals emitted by each of said plurality of sensors and said concentrations of sterilant vapor for each of said plurality of sensors.

3. The method of claim 2, wherein temperature of air ranges from about -10 °C to about 85 °C.

- 5 4. The method of claim 2, wherein relative humidity of air ranges from about 10 to about 100%.
 - 5. The method of claim 2, wherein flow rate of air ranges from about 1 to about 200 cubic meters per hour.
 - 6. The method of claim 2, wherein concentration of sterilant ranges from 100 to 25000 ppm.
 - 7. The method of claim 2, wherein said sterilant is hydrogen peroxide.
 - 8. The method of claim 1, wherein said sensor calibration data is generated by a method comprising the steps of:
 - (a) providing a sensor;

10

15

- 20 (b) subjecting said sensor to at least two quantities of air, each of said at least two quantities of air having (1) a known quality and (2) a known concentration of sterilant vapor, said sterilant vapor having a known physical condition;
 - (c) measuring signals emitted by said sensor, each of said signals corresponding to a concentration of sterilant vapor; and
 - (d) establishing a mathematical relationship between said signals emitted and said concentrations of sterilant vapor for said sensor.
- 9. The method of claim 8, wherein temperature of air ranges from about 30 -10 °C to about 85 °C.
 - 10. The method of claim 8, wherein relative humidity of air ranges from about 10 to about 100%.

11. The method of claim 8, wherein flow rate of air ranges from about 1 to about 200 cubic meters per hour.

- 12. The method of claim 8, wherein concentration of sterilant ranges from 100 to 25000 ppm.
 - 13. The method of claim 8, wherein said sterilant is hydrogen peroxide.
- The method of claim 1, wherein said sensor calibration data is
 normalized to compensate for the difference between said measurable parameter for said reference calibration data and said measurable parameter for said sensor calibration data for said sensor by a method comprising the steps of:
 - (a) selecting a concentration of sterilant vapor;

15

20

- (b) determining the value of said measurable parameter at which said concentration of said sterilant vapor obtained from said sensor calibration data equals said concentration of said sterilant vapor obtained from said reference calibration data;
- (c) adjusting the values measured by said individual sensor a sufficient amount to compensate for the deviation between said reference calibration data and said sensor calibration data.
- 15. A method for calibrating a sensor that is used for measuring the quantity of a sterilant in a system for delivering said sterilant, said method comprising the steps of:
 - (a) preparing a reference calibration curve, said reference calibration curve having a slope and an intercept;
 - (b) preparing a sensor calibration curve, said sensor calibration curve having a slope and an intercept;
 - (c) normalizing said sensor calibration curve to compensate for (1) the difference between said slope of said reference calibration curve and said slope of said sensor calibration curve and (2) the difference between said

intercept of said reference calibration curve and said intercept of said sensor calibration curve.

- 16. The method of claim 15, wherein said reference calibration curve is prepared by a method comprising the steps of:
 - (a) providing a plurality of sensors;

10

15

20

- (b) subjecting each of said plurality of sensors to at least two quantities of air, each of said at least two quantities of air having (i) a known quality and (ii) a known concentration of sterilant vapor, said sterilant vapor having a known physical condition;
- (c) measuring said signals emitted by each of said plurality of sensors, each of said signals being proportional to said concentration of sterilant vapor;
- (d) establishing a linear mathematical relationship between said signals emitted by each of said plurality of sensors and said concentrations of sterilant vapor for each of said plurality of sensors; and
- (e) establishing said reference curve by means of a statistical analysis of said signals emitted by each of said plurality of sensors and said concentrations of sterilant vapor for each of said plurality of sensors.
- 17. The method of claim 16, wherein temperature of air ranges from about 10 °C to about 85 °C.
- The method of claim 16, wherein relative humidity of air ranges from about 10 to about 100%.
 - 19. The method of claim 16, wherein flow rate of air ranges from about 1 to about 200 cubic meters per hour.
 - 20. The method of claim 16, wherein concentration of sterilant ranges from 100 to 25000 ppm.
 - 21. The method of claim 16, wherein said sterilant is hydrogen peroxide.

22. The method of claim 15, wherein said sensor calibration curve is prepared by a method comprising the steps of:

(a) providing an individual sensor;

5

10

15

- (b) subjecting said sensor to at least two quantities of air, each of said at least two quantities of air having (1) a known quality and (2) a known concentration of sterilant vapor, said sterilant vapor having a known physical condition;
- (c) measuring said signals emitted by said sensor, each of said signals being proportional to said concentration of sterilant vapor; and
- (d) establishing a linear mathematical relationship between said signals emitted by said individual sensor and said concentrations of sterilant vapor for said individual sensor.

23. The method of claim 22, wherein temperature of air ranges from about - 10 °C to about 85 °C.

- 24. The method of claim 22, wherein relative humidity of air ranges from about 10 to about 100%.
 - 25. The method of claim 22, wherein flow rate of air ranges from about 1 to about 200 cubic meters per hour.
- 25 26. The method of claim 22, wherein concentration of sterilant ranges from 100 to 25000 ppm.
 - 27. The method of claim 22, wherein said sterilant is hydrogen peroxide.
- 28. The method of claim 15, wherein said sensor calibration curve is normalized to compensate for (1) the difference between said slope of said reference calibration curve and said slope of said sensor calibration curve and (2) the difference between said intercept of said reference calibration curve and said intercept of said sensor calibration curve by a method comprising the steps of:

(a) determining said intercept of said reference calibration curve;

- (c) determining said intercept of said sensor calibration curve;
- (d) determining said slope of said reference calibration curve;
- (e) determining said slope of said sensor calibration curve;

5

10

15


- (f) adjusting said sensor calibration curve, if necessary, in order to compensate for said difference between said intercept of said reference calibration curve and said intercept of said sensor calibration curve;
- (g) adjusting said sensor calibration curve, if necessary, in order to compensate for said difference between said slope of said reference calibration curve and said slope of said sensor calibration curve.
- 29. A method for calibrating a portable unit for measuring the concentration of a sterilant, said method comprising the steps of:
 - (a) passing a sterilant vapor through a vessel, said vessel being submerged in a water bath, said portable unit being installed within said vessel, said portable unit responsive to the concentration of sterilant, the temperature, and the relative humidity for a test run; and
 - (b) determining the concentration of said sterilant by a titration method.
 - 30. The method of claim 29, wherein said sterilant is hydrogen peroxide.
- 31. A method for verifying that a sensor for measuring the quantity of sterilant in a system for delivering said sterilant has been calibrated accurately, wherein said sterilant is delivered in a stream of flowing air, said method comprising the steps of:
- (a) performing at least two static test runs for said sensor to obtain readings of concentration of sterilant vapor;
 - (b) performing at least two dynamic test runs for said sensor to obtain readings of concentration of sterilant vapor;
 - (c) comparing the results of said at least two static test runs;
 - (d) comparing the results of said at least two dynamic test runs;

(e) concluding that (1) the flow rate of air has changed or (2) the rate of addition of sterilant to the system has changed or (3) said sensor is unreliable if (i) if said readings of concentration of sterilant vapor determined in said static test runs remain substantially constant, but said readings of concentration of sterilant vapor determined in said dynamic test runs vary substantially or (ii) if said readings of concentration of sterilant vapor determined in said static test runs vary substantially, but said readings of concentration of sterilant vapor determined in said dynamic test runs remain substantially constant.


5


10

32. The method of claim 31, wherein said sterilant is hydrogen peroxide.

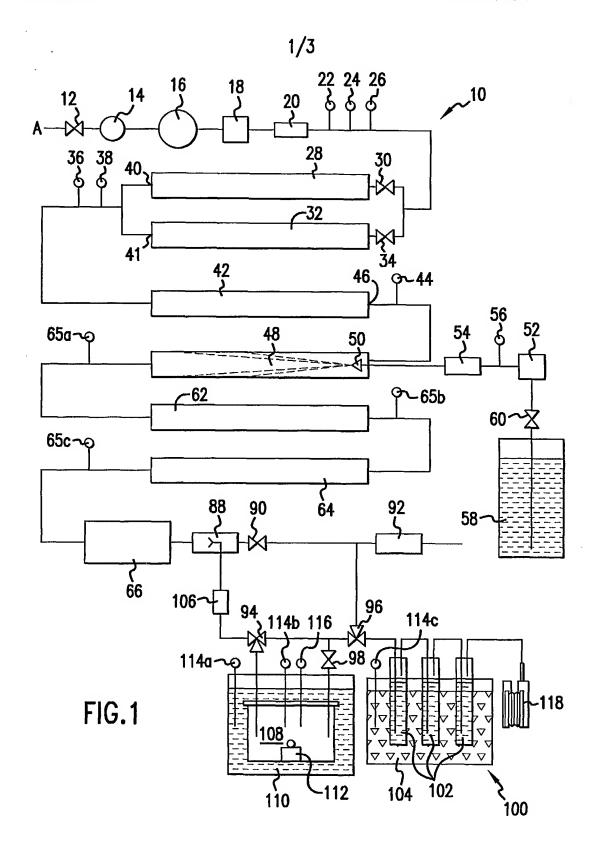

.

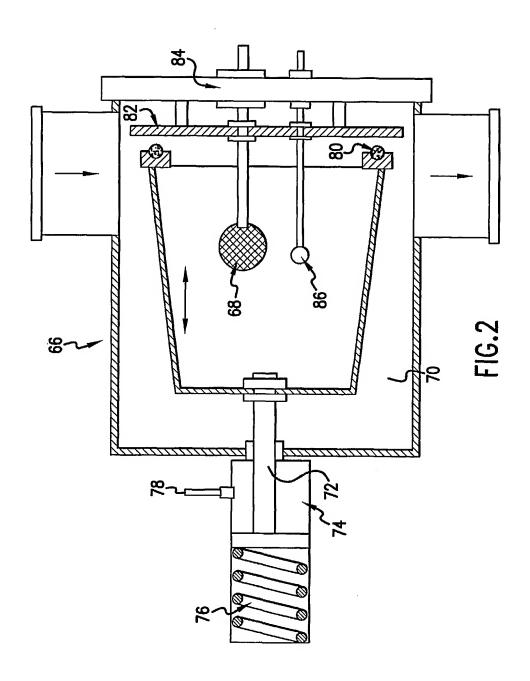

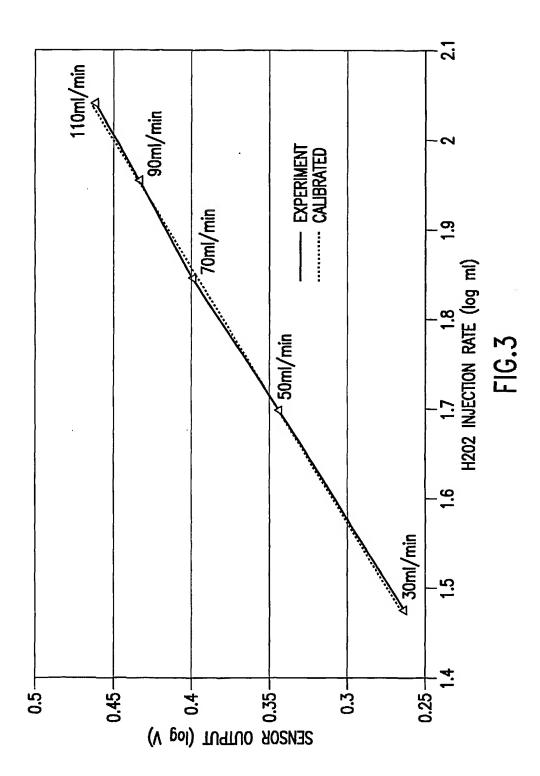

INTERNATIONAL SEARCH REPORT

PCT/DK 03/00728

A. CLASSI IPC 7	FICATION OF SUBJECT MATTER G01N27/00 G01N33/00			
0.0004500.00	Later referred Detect Classification (IDO) and the United Section (IDO)			
	o International Patent Classification (IPC) or to both national classification (IPC) or to both national classification (IPC)	ation and IPC		
	ocumentation searched (classification system followed by classification	on symbols)	(F) 0.	1
IPC 7	G01N			المال المال المال
Documental	ion searched other than minimum documentation to the extent that s	uch documents are included in the fields sea	-1, 2	
Electronic d	ata base consulted during the international search (name of data ba	se and, where practical, search terms used)	Copol. Co	Pennanco
EPO-In	ternal, WPI Data, PAJ, INSPEC			
C. DOCUMI	ENTS CONSIDERED TO BE RELEVANT			
Category °	Citation of document, with indication, where appropriate, of the rela	evant passages	Relevant to claim No.	
А	US 6 037 178 A (LEINER MARCO JEA ET AL) 14 March 2000 (2000-03-14 column 3, line 45 -column 4, line abstract)	1-23	
А	WO 02 <u>066973</u> A (ABBOTT LAB) 29 August 2002 (2002-08-29) abstract		1-23	
A	US 4 163 734 A (HOLBEK CARL C E 7 August 1979 (1979-08-07) abstract	T AL)	1-23	
A	US 4 786 394 A (BURGESS BRUCE M 22 November 1988 (1988-11-22) abstract; figure 1	ET AL)	1-23	
		-/		
X Furth	er documents are listed in the continuation of box C.	X Patent family members are listed in	n annex.	
	regories of cited documents :			
"A" dacume conside	art defining the general state of the art which is not ered to be of particular relevance	"T" later document published after the Inter or priority date and not in conflict with it cited to understand the principle or the invention	ory underlying the	
filing da	ate	"X" document of particular relevance; the cl cannot be considered novel or cannot involve an inventive step when the doc	be considered to	
which i	nt which may throw doubts on priority daim(s) or s cited to establish the publication date of another or other special reason (as specified)	"Y" document of particular relevance; the ci-	aimed Invention	
	nt referring to an oral disclosure, use, exhibition or	document is combined with one or mor ments, such combination being obviou	re other such docu-	
"P" docume	int published prior to the international filing date but an the priority date claimed	in the art. "&" document member of the same patent f		
Date of the a	ctual completion of the international search	Date of mailing of the international sear	rch report	1
29	January 2004	0 9. 02. 04		
Name and m	ailing address of the ISA	Authorized officer	-	
•	European Patent Office, P.B. 5818 Patentlaan 2 NL - 2280 HX Rijswijk Tel. (+31-70) 340-2040, Tx. 31 651 epo nl, Fax: (+31-70) 340-3016	GORDANA NINKOVIC/	JAA	

INTERNATIONAL SEARCH REPORT

Interns oplication No
PCT/FK 03/00728

		PCT/bk 03/00728			
	(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT				
Category °	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.			
4	US 4 469 792 A (OWEN WELLS P ET AL) 4 September 1984 (1984-09-04) abstract	1-23			
	•				
į					
ļ					

INTERNATIONAL SEARCH REPORT

mation on patent family members

PCT/TK 03/00728

Patent document cited in search report	Publication date	Patent far member		Publication date
US 6037178 A	14-03-2000	AT 121 DE 59604 EP 0754 JP 2928	1944 A1	26-02-2001 15-06-2000 17-02-2000 22-01-1997 03-08-1999 07-02-1997
WO 02066973 A	29-08-2002		7711 A1 2238 A2 3616 A 5973 A2	24-10-2002 29-08-2002 19-11-2003 15-10-2003 29-08-2002 08-05-2003
US 4163734 A	07-08-1979	US 4116	1975 A 5336 A L108 A	28-02-1977 26-09-1978 24-04-1979
US 4786394 A	22-11-1988	CA 1246 CA 1313 DE 3683 EP 0226 JP 7051	1726 T 5674 A1 3486 C2 3471 D1 5593 A1 1127 B 5590 A1	15-02-1992 13-12-1988 09-02-1993 27-02-1992 01-07-1987 05-06-1995 25-09-1986
US 4469792 A	04-09-1984	US 4401	1652 A	30-08-1983